

Features:

Jumbo LCD

Accuracy:  $\pm 0.05^{\circ}\text{C}$  ( $\pm 0.1^{\circ}\text{F}$ )

Fever alarm

Last reading memory

Rigid tip

Low battery indicator

Automatic Off

Designed for oral use




User Instructions:

Step1. Press the On/Off Button next to LCD display. A tone will sound as the screen shows , followed by last memory temperature.

After showing the self-test temperature, the thermometer is now in the testing mode.

\*Note: If the measured temperature is less than  $32^{\circ}\text{C}$  or  $90^{\circ}\text{F}$ , the LCD will display Lo, if the measured temperature is higher than  $42.9^{\circ}\text{C}$  or  $109.9^{\circ}\text{F}$ , the LCD will display Hi.

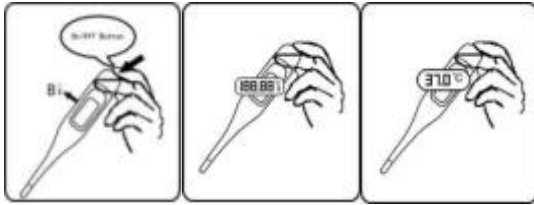


Step2. Position thermometer in desired location (mouth)

a) Oral Use: Place thermometer under tongue as indicated by "v" position. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air. Normal temperature between 35.7°C and 37.3°C (96.3°F and 99.1°F)

Step3.

The degree sign flashes throughout the testing process. When flashing stops an alarm will beep for approximately 10 seconds. The measured reading will appear on the LCD simultaneously. The minimum measurement time until the signaling tone (beep) must be maintained without exception. The measurement continues even after the buzzer notification. The displayed temperature will not change when thermometer is removed from its testing position.

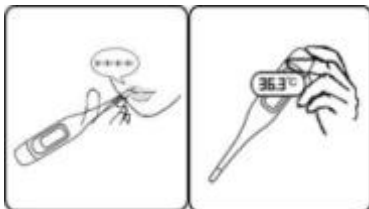

Step4.

To prolong battery life, press the On/Off Button to turn unit off after testing is complete. If no action is taken, the unit will automatically shut off after around 10 minutes. Store the thermometer in its protective case.

